

# PHYSICS PRACTICAL SHEETS

.....KV..... Campus

Date 11 June 2026

Class : CE - 1/11

Roll No.: 43

Shift: Morning

Object of the Experiment (Block Letter)

*Rausulfergi*

Experiment No.: 2

Group : S

Sub: Physics lab 102

Set: 1

## MEASUREMENT OF MOMENT OF INERTIA OF FLYWHEEL

Apparatus required:

- 1) Flywheel
- 2) Stopwatch
- 3) Meter Scale
- 4) Thread
- 5) weight hanger
- 6) known weights
- 7) Marker
- 8) Vernier Calipers

Theory:

When a mass  $m$ , attached to the axle of a flywheel via a thread, falls through a height  $h$ , its potential energy is transferred into:

- a) Translational kinetic energy of the mass
- b) Rotational kinetic energy of the flywheel
- c) Work done against friction

Using the principles of energy conservation and accounting for the work done against friction, the moment of inertia ( $I$ ) of the flywheel can be determined using the following formula:

$$I = \frac{2mgh - mr^2\omega^2}{\omega^2 \left(1 + \frac{n}{n_1}\right)} \quad \text{--- (1)}$$



where,

$r$  is the radius of the axle

$n$  is number of turns of the thread wound around the axle

$n_1$  is the number of revolutions made by the flywheel before coming to rest after the string is detached.

$\omega$  is the angular velocity of the flywheel

$$\omega = \frac{2 \times 2\pi n_1}{t} \quad \text{--- (2)}$$

$h$  is the vertical distance through which the mass  $m$  falls as the string unwinds from axle of the flywheel

Observations:

Vernier constant =  $0.1 \text{ mm} = 0.01 \text{ cm}$

Diameter of the axle,  $d =$  (i)  $2.14 \text{ cm}$  (ii)  $2.16 \text{ cm}$  (iii)  $2.15 \text{ cm}$

Mean diameter ( $d$ ) =  $2.15 \text{ cm}$

Radius of the axle,  $r = 1.075 \text{ cm}$

Circumference of the wheel =  $58 \text{ cm}$

For mass: (i)  $m = 200 \text{ g}$ , Height  $h = 71 \text{ cm}$

No. of turns wound on axle,  $n = 10$

(ii)  $m = 250 \text{ g}$ , height  $h = 70.5 \text{ cm}$

No. of turns wound on axle,  $n = 10$

(iii)  $m = 300 \text{ g}$ , height  $h = 71 \text{ cm}$

No. of turns wound on axle,  $n = 10$



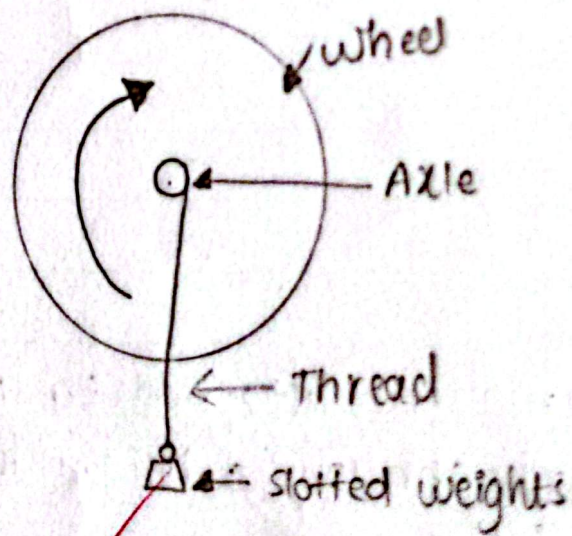
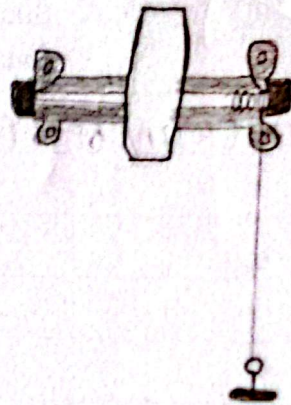


Fig 1: Flywheel mounted on the wall



Table 1: Determination of the moment of inertia of the flywheel

| Mass<br>m | S<br>N | No of complete<br>revs (made by<br>the wheel) x | Fraction<br>of revolution | Total<br>No. of<br>revs (x+y) | Mean<br>$n_1$ | Time | Mean<br>time<br>t | Angular<br>velocity<br>( $\omega$ )<br>$\omega = \frac{2\pi n_1}{t}$ | Moment of<br>inertia<br>$I = 2mgh - mr^2\omega^2$<br>$\omega^2 \left(1 + \frac{r^2}{n_1^2}\right)$ |
|-----------|--------|-------------------------------------------------|---------------------------|-------------------------------|---------------|------|-------------------|----------------------------------------------------------------------|----------------------------------------------------------------------------------------------------|
| 200<br>gm | 1      | 14                                              | 0.71                      | 14.71                         | 14.74         | 27   |                   |                                                                      |                                                                                                    |
|           | 2      | 14                                              | 0.69                      | 14.69                         |               | 28   | 26.67             | 6.94                                                                 | $3.44 \times 10^5$                                                                                 |
|           | 3      | 14                                              | 0.82                      | 14.82                         |               | 25   |                   |                                                                      |                                                                                                    |
| 250<br>gm | 1      | 20                                              | 0.09                      | 20.09                         | 19.58         | 30   | 30.3              | 8.12                                                                 | $3.46 \times 10^5$                                                                                 |
|           | 2      | 19                                              | 0.01                      | 19.01                         |               | 28   |                   |                                                                      |                                                                                                    |
|           | 3      | 19                                              | 0.56                      | 19.56                         |               | 33   |                   |                                                                      |                                                                                                    |
| 300<br>gm | 1      | 22                                              | 0.76                      | 22.76                         | 23.12         | 30   | 31                | 9.37                                                                 | $3.31 \times 10^5$                                                                                 |
|           | 2      | 23                                              | 0.77                      | 23.77                         |               | 32   |                   |                                                                      |                                                                                                    |
|           | 3      | 22                                              | 0.83                      | 22.83                         |               | 31   |                   |                                                                      |                                                                                                    |

Mean moment of inertia,  $I = 3.4 \times 10^5 \text{ gm cm}^2$

Result:

The moment of Inertia (I) of the flywheel is found to be  $3.4 \times 10^5 \text{ gm cm}^2$

### Precautions:

- 1) Minimize friction in the flywheel and its axle to ensure accurate measurements.
- 2) Check that the string is free of kinks or twists to avoid uneven rotation.
- 3) Use a thin, strong string & wind it evenly around the axle for consistent results.
- 4) Start the stopwatch precisely when the mass detaches from the axle to record accurate timing.

Conclusion: The moment of inertia of the flywheel is determined from the experiment.

*Pranav Singh*  
11th June 2026.